

## Total RNA Purification Kit Dx

**REF** Dx17200



**IVD**

## Product Insert

 PIDx17200-1

### Intended Use

Norgen's Total RNA Purification Kit Dx provides a rapid method for the isolation and purification of total RNA from tissue samples, blood, plasma, serum, bacteria, yeast, fungi and viruses for subsequent *in vitro* diagnostic use. The kit purifies all sizes of RNA, from large mRNA and ribosomal RNA down to microRNA (miRNA) and small interfering RNA (siRNA). The RNA is preferentially purified from other cellular components such as proteins, without the use of phenol or chloroform.

This kit is designed to be used with any downstream application employing enzymatic amplification or other enzymatic modifications of RNA followed by signal detection or amplification. Any diagnostic results generated using the RNA isolated with Norgen's Total RNA Purification Kit Dx in conjunction with an *in vitro* diagnostic assay should be interpreted with regard to other clinical or laboratory findings.

To minimize irregularities in diagnostic results, suitable controls for downstream applications should be used.










Norgen's Total RNA Purification Kit Dx is intended for use by professional users such as technicians, physicians and biologists experienced and trained in molecular biological techniques.

Norgen's Total RNA Purification Kit Dx does not provide a diagnostic result. It is the sole responsibility of the user to use and validate the kit in conjunction with a downstream *in vitro* diagnostic assay.

### Kit Components

Component	Product # Dx17200 (50 preps)
Lysis Solution	40 mL
Wash Solution	22 mL
Elution Solution	6 mL
Mini Spin Columns	50
Collection Tubes	50
Elution tubes (1.7 mL)	50
Product Insert	1

### Label Legend

								
Do not reuse	Use by	Batch Code	Catalogue Number	Contains sufficient for <n> tests	Manufacturer	In Vitro Diagnostic Medical Device	Consult instructions for use	Temperature limitation

### Advantages

- CE-IVD marked in accordance with EU Directive 98/79/EC
- Fits into *in vitro* diagnostic workflows
- Fast and easy processing using rapid spin-column format
- Isolate total RNA, from large rRNA down to microRNA (miRNA)
- No phenol or chloroform extractions
- Isolate high quality total RNA from a variety of sources
- RNA can be isolated and detected from as little as a single animal cell

### Specifications

Kit Specifications	
Maximum Column Binding Capacity	50 µg
Maximum Column Loading Volume	650 µL
Size of RNA Purified	All sizes, including small RNA (<200 nt)
Maximum Amount of Starting Material:	
Tissues	10 mg (for most tissues)
Blood	100 µL
Plasma/Serum	200 µL
Bacteria	1 x 10 <sup>9</sup> cells
Yeast	1 x 10 <sup>8</sup> cells
Fungi	50 mg
Time to Complete 10 Purifications	20 minutes
Average Yields	
HeLa Cells (1 x 10 <sup>6</sup> cells)	15 µg
<i>E. coli</i> (1 x 10 <sup>9</sup> cells)	50 µg

### Storage Conditions and Product Stability

All solutions should be kept tightly sealed and stored at room temperature. All solutions and plastics can be used until the expiration date specified on their labels.

### Precautions

Ensure that a suitable lab coat, disposable gloves and protective goggles are worn when working with chemicals. For more information, please consult the appropriate Material Safety Data Sheets (MSDSs). These are available as convenient PDF files online at [www.norgenbiotek.com](http://www.norgenbiotek.com).

The **Lysis Solution** contains guanidinium salts, and should be handled with care. Guanidinium salts form highly reactive compounds when combined with bleach, thus care must be taken to properly dispose of any of these solutions

Blood of all human and animal subjects is considered potentially infectious. All necessary precautions recommended by the appropriate authorities in the country of use should be taken when working with whole blood.

## Customer-Supplied Reagents and Equipment

You must have the following in order to use the Total RNA Purification Kit Dx:

### *For All Protocols*

- Benchtop microcentrifuge
- 96 - 100% ethanol
- $\beta$ -mercaptoethanol (optional)

### *For Tissue Protocol*

- Liquid nitrogen
- Mortar and pestle
- 70% ethanol

### *For Nasal or Throat Swabs*

- Sterile, single-use cotton swabs

### *For Bacterial Protocol*

- Lysozyme-containing TE Buffer:
  - For Gram-negative bacteria, 1 mg/mL lysozyme in TE Buffer
  - For Gram-positive bacteria, 3 mg/mL lysozyme in TE Buffer

### *For Yeast Protocol*

- Resuspension Buffer with Lyticase:
  - 50 mM Tris pH 7.5
  - 10 mM EDTA
  - 1 M Sorbitol
  - 1 unit/ $\mu$ L Lyticase

### *For Fungi Protocol*

- Liquid nitrogen
- Mortar and pestle
- 70% ethanol

### *For Plasma/Serum Protocol*

- MS2 RNA (0.8  $\mu$ g/ $\mu$ l). (Roche, Cat. No. 10165948001)

## Working with RNA

RNases are very stable and robust enzymes that degrade RNA. Autoclaving solutions and glassware is not always sufficient to actively remove these enzymes. The first step when preparing to work with RNA is to create an RNase-free environment. The following precautions are recommended as your best defense against these enzymes.

- The RNA area should be located away from microbiological work stations
- Clean, disposable gloves should be worn at all times when handling reagents, samples, pipettes, disposable tubes, etc. It is recommended that gloves are changed frequently to avoid contamination
- There should be designated solutions, tips, tubes, lab coats, pipettes, etc. for RNA only
- All RNA solutions should be prepared using at least 0.05% DEPC-treated autoclaved water or molecular biology grade nuclease-free water
- Clean all surfaces with commercially available RNase decontamination solutions
- When working with purified RNA samples, ensure that they remain on ice during downstream applications

## Procedures


All centrifugation steps are carried out in a benchtop microcentrifuge. Various speeds are required for different steps, so please check your microcentrifuge specifications to ensure that it is capable of the proper speeds. All centrifugation steps are performed at room temperature. The correct rpm can be calculated using the formula:

$$RPM = \sqrt{\frac{RCF}{(1.118 \times 10^{-5}) (r)}}$$

where  $RCF$  = required gravitational acceleration (relative centrifugal force in units of  $g$ );  $r$  = radius of the rotor in  $cm$ ; and  $RPM$  = the number of revolutions per minute required to achieve the necessary  $g$ -force.

## Flow Chart

Procedure for Purifying Total RNA using Norgen's Total RNA Purification Kit Dx

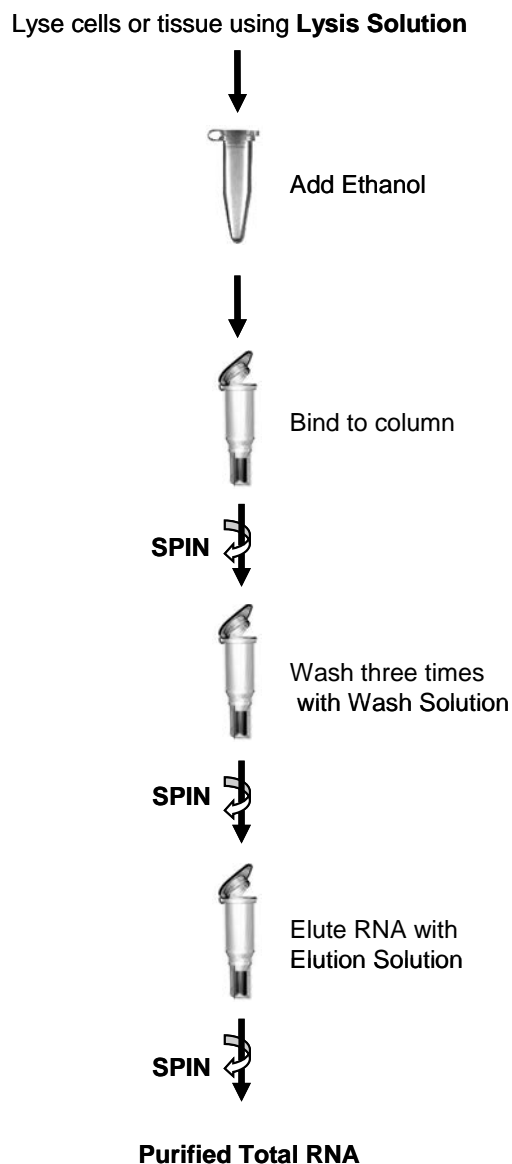

## Section 1. Preparation of Lysate From Various Cell Types

### Notes Prior to Use

- *The steps for preparing the lysate are different depending on the starting material (Step 1). However, the subsequent steps are the same in all cases (Steps 2 – 6), with the exception of the protocol for plasma/serum. A separate protocol for the isolation of total RNA from plasma/serum samples is located in Appendix B.*
- Please ensure that the correct procedure for preparing the lysate from your starting material is followed, as indicated in the table below:

Sample Type	Lysate Preparation Page #
Tissue	5
Blood	6
Plasma/Serum	11
Nasal / Throat Swabs	7
Bacteria	7
Yeast	8
Fungi	8
Viruses	9

- All centrifugation steps are carried out in a benchtop microcentrifuge at 14,000 x g (~ 14,000 RPM) except where noted. All centrifugation steps are performed at room temperature.
- A variable speed centrifuge should be used for maximum kit performance. If a variable speed centrifuge is not available a fixed speed centrifuge can be used, however reduced yields may be observed.
- Ensure that all solutions are at room temperature prior to use.
- Prepare a working concentration of the **Wash Solution** by adding 50 mL of 96-100% ethanol (provided by the user) to the supplied bottle containing the concentrated **Wash Solution**. This will give a final volume of 72 mL. The label on the bottle has a box that may be checked to indicate that the ethanol has been added.
- **Optional:** The use of  $\beta$ -mercaptoethanol in lysis is highly recommended for most tissues, particularly those known to have a high RNase content (ex: pancreas), as well as for nasal and throat swabs. It is also recommended for users who wish to isolate RNA for sensitive downstream applications. Add 10  $\mu$ L of  $\beta$ -mercaptoethanol (provided by the user) to each 1 mL of Lysis Solution required.  $\beta$ -mercaptoethanol is toxic and should be dispensed in a fume hood. Alternatively, the lysis solution can be used as provided.
- It is important to work quickly during this procedure.

### 1A. Lysate Preparation from Tissues

#### Notes Prior to Use

- Norgen's Total RNA Purification Kit is designed for isolating RNA from small amount of tissue sample (up to 10 mg in most cases). If a larger amount of starting material is desired, Norgen's Animal Tissue RNA Purification Kit (Cat.# 25700) should be used.
- RNA in tissues is not protected after harvesting until it is disrupted and homogenized. Thus it is important that the procedure is carried out as quickly as possible, particularly the Cell Lysate Preparation step.
- Fresh or frozen tissues may be used for the procedure. Tissues should be flash-frozen in liquid nitrogen and transferred immediately to a -70°C freezer for long-term storage. Tissues may be stored at -70°C for several months. When isolating total RNA from

- frozen tissues ensure that the tissue does not thaw during weighing or prior to grinding with the mortar and pestle.
- Tissues stored in RNA stabilization reagents such as RNA*later*<sup>®</sup> are compatible with this isolation procedure. Prior to isolation, carefully remove the tissue from the storage reagent using forceps, and dry excessive liquid.
  - The maximum recommended input of tissue varies depending on the type of tissue being used. Please refer to Table 1 below as a guideline for maximum tissue input amounts. If your tissue of interest is not included in the table below we recommend starting with an input of no more than 10 mg.

**Table 1. Recommended Maximum Input Amounts of Different Tissues**

Tissue	Maximum Input Amount
Brain	25 mg
Heart	5 mg
Kidney	10 mg
Liver	10 mg
Lung	10 mg
Spleen	10 mg

#### 1A. Cell Lysate Preparation from Tissues

- Excise the tissue sample from the animal.
- Determine the amount of tissue by weighing. Please refer to Table 1 for the recommended maximum input amounts of different tissues. For tissues not included in the table, we recommend starting with an input of no more than 10 mg.
- Transfer the tissue into a mortar that contains an appropriate amount of liquid nitrogen to cover the sample. Grind the tissue thoroughly using a pestle.
- Allow the liquid nitrogen to evaporate, without allowing the tissue to thaw.
- Add 600  $\mu$ L of **Lysis Solution** to the tissue sample and continue to grind until the sample has been homogenized. Homogenize by passing the lysate 5-10 times through a 25 gauge needle attached to a syringe.
- Using a pipette, transfer the lysate into an RNase-free microcentrifuge tube (not provided).
- Spin the lysate for 2 minutes to pellet any cell debris. Transfer the supernatant to another RNase-free microcentrifuge tube. Note the volume of the supernatant/lysate.
- Add an equal volume of 70% ethanol (provided by the user) to the lysate volume collected (100  $\mu$ L of ethanol is added to every 100  $\mu$ L of lysate). Vortex to mix. **Proceed to Step 2 (page 9).**

#### 1B. Lysate Preparation from Blood

##### Notes Prior to Use

- This procedure is for the isolation of RNA from whole blood. **For the isolation of RNA from plasma or serum samples, please see Appendix B.**
- Blood of all human and animal subjects is considered potentially infectious. All necessary precautions recommended by the appropriate authorities in the country of use should be taken when working with whole blood.
- It is recommended that no more than 100  $\mu$ L of blood be used in order to prevent clogging of the column.

- We recommend the use of this kit to isolate RNA from non-coagulating fresh blood using EDTA as the anti-coagulant.
- It is important to work quickly during this procedure.

### 1B. Cell Lysate Preparation from Blood

- a. Transfer up to 100  $\mu\text{L}$  of non-coagulating blood to an RNase-free microcentrifuge tube (not provided).
- b. Add 350  $\mu\text{L}$  of **Lysis Solution** to the blood. Lyse cells by vortexing for 15 seconds. Ensure that mixture becomes transparent before proceeding to the next step.
- c. Add 200  $\mu\text{L}$  of 95 – 100% ethanol (provided by the user) to the lysate. Mix by vortexing for 10 seconds. **Proceed to Step 2 (page 9).**

### 1C. Lysate Preparation from Nasal or Throat Swabs

#### Notes Prior to Use

- Body fluid of all human and animal subjects is considered potentially infectious. All necessary precautions recommended by the appropriate authorities in the country of use should be taken when working with these samples.
- It is important to work quickly during this procedure.

### 1C. Cell Lysate Preparation from Nasal or Throat Swabs

- a. Add 600  $\mu\text{L}$  of **Lysis Solution** to an RNase-free microcentrifuge tube (not provided).
- b. Gently brush a sterile, single-use cotton swab inside the nose or mouth of the subject.
- c. Using sterile techniques, cut the cotton tip where the nasal or throat cells were collected and place into the microcentrifuge tube containing the **Lysis Solution**. Close the tube. Vortex gently and incubate for 5 minutes at room temperature.
- d. Using a pipette, transfer the lysate into another RNase-free microcentrifuge tube (not provided). Note the volume of the lysate.
- e. Add an equal volume of 70% ethanol (provided by the user) to the lysate volume collected (100  $\mu\text{L}$  of ethanol is added to every 100  $\mu\text{L}$  of lysate). Vortex to mix. **Proceed to Step 2 (page 9).**

### 1D. Lysate Preparation from Bacteria

#### Notes Prior to Use

- Prepare the appropriate lysozyme-containing TE Buffer as indicated in Table 2. This solution should be prepared with sterile, RNase-free TE Buffer, and kept on ice until needed. These reagents are to be provided by the user.
- It is recommended that no more than  $10^9$  bacterial cells be used in this procedure. Bacterial growth can be measured using a spectrophotometer. As a general rule, an *E. coli* culture containing  $1 \times 10^9$  cells/mL has an  $\text{OD}_{600}$  of 1.0.
- For RNA isolation, bacteria should be harvested in log-phase growth.
- Bacterial pellets can be stored at  $-70^\circ\text{C}$  for later use, or used directly in this procedure.
- Frozen bacterial pellets should not be thawed prior to beginning the protocol. Add the Lysozyme-containing TE Buffer directly to the frozen bacterial pellet (**Step 1Ec**).

### 1D. Cell Lysate Preparation from Bacteria

- a. Pellet bacteria by centrifuging at  $14,000 \times g$  (~14,000 RPM) for 1 minute.
- b. Decant supernatant, and carefully remove any remaining media by aspiration.
- c. Resuspend the bacteria thoroughly in 100  $\mu\text{L}$  of the appropriate lysozyme-containing TE buffer (see Table 1) by vortexing. Incubate at room temperature for the time indicated in Table 1.

- d. Add 300  $\mu$ L of Lysis Solution and vortex vigorously for at least 10 seconds.
- e. Add 200  $\mu$ L of 95 – 100% ethanol (provided by the user) to the lysate. Mix by vortexing for 10 seconds. **Proceed to Step 2 (page 9).**

**Table 2: Incubation Time for Different Bacterial Strains**

Bacteria Type	Lysozyme Concentration in TE Buffer	Incubation Time
Gram-negative	1 mg/mL	5 min
Gram-positive	3 mg/mL	10 min

## 1E. Lysate Preparation from Yeast

### Notes Prior to Use

- Prepare the appropriate amount of Lyticase-containing Resuspension Buffer, considering that 100  $\mu$ L of buffer is required for each preparation. The Resuspension Buffer should have the following composition: 50 mM Tris, pH 7.5, 10 mM EDTA, 1M Sorbital, 0.1%  $\beta$ -mercaptoethanol and 1 unit/ $\mu$ L Lyticase. This solution should be prepared with sterile, RNase-free reagents, and kept on ice until needed. These reagents are to be provided by the user.
- It is recommended that no more than  $10^7$  yeast cells or 1 mL of culture be used for this procedure.
- For RNA isolation, yeast should be harvested in log-phase growth.
- Yeast can be stored at  $-70^{\circ}\text{C}$  for later use, or used directly in this procedure.
- Frozen yeast pellets should not be thawed prior to beginning the protocol. Add the Lyticase-containing Resuspension Buffer directly to the frozen yeast pellet (**Step 1Fc**).

### 1E. Cell Lysate Preparation

- a. Pellet yeast by centrifuging at  $14,000 \times g$  (~14,000 RPM) for 1 minute.
- b. Decant supernatant, and carefully remove any remaining media by aspiration.
- c. Resuspend the yeast thoroughly in 100  $\mu$ L of Lyticase-containing Resuspension Buffer by vortexing. Incubate at  $37^{\circ}\text{C}$  for 10 minutes.
- d. Add 300  $\mu$ L of Lysis Solution and vortex vigorously for at least 10 seconds.
- e. Add 200  $\mu$ L of 95 – 100% ethanol (provided by the user) to the lysate. Mix by vortexing for 10 seconds. **Proceed to Step 2 (page 9).**

## 1F. Lysate Preparation from Fungi

### Notes Prior to Use

- Fresh or frozen fungi may be used for this procedure. Fungal tissues should be flash-frozen in liquid nitrogen and transferred immediately to a  $-70^{\circ}\text{C}$  freezer for long-term storage. Fungi may be stored at  $-70^{\circ}\text{C}$  for several months. Do not allow frozen tissues to thaw prior to grinding with the mortar and pestle in order to ensure that the integrity of the RNA is not compromised.
- It is recommended that no more than 50 mg of fungi be used for this procedure in order to prevent clogging of the column.

### 1F. Cell Lysate Preparation from Fungi

- a. Determine the amount of fungi by weighing. It is recommended that no more than 50 mg of fungi be used for the protocol.



- b. Transfer the fungus into a mortar that contains an appropriate amount of liquid nitrogen to cover the sample. Grind the fungus thoroughly using a pestle.

**Note:** At this stage the ground fungus may be stored at  $-70^{\circ}\text{C}$ , such that the RNA purification can be performed at a later time.

- c. Allow the liquid nitrogen to evaporate, without allowing the tissue to thaw.
- d. Add 600  $\mu\text{L}$  of **Lysis Solution** to the tissue sample and continue to grind until the sample has been homogenized.
- e. Using a pipette, transfer the lysate into an RNase-free microcentrifuge tube (not provided).
- f. Spin the lysate for 2 minutes to pellet any cell debris. Transfer the supernatant to another RNase-free microcentrifuge tube. Note the volume of the supernatant/lysate.
- g. Add an equal volume of 70% ethanol (provided by the user) that is equivalent to the lysate volume collected (100  $\mu\text{L}$  of ethanol is added to every 100  $\mu\text{L}$  of lysate). Vortex to mix. **Proceed to Step 2 (page 9).**

## 1G. Lysate Preparation from Viruses

### *Notes Prior to Use*

- For the isolation of integrated viral RNA, follow Section **1A** if the starting material is cell culture, follow Section **1B** if the starting material is tissue, follow Section **1C** if the starting material is blood, or follow Section **1H** if the starting material is a nasal or throat swab.
- For the isolation of RNA from free viral particles, follow the procedure below.
- It is recommended that no more than 100  $\mu\text{L}$  of viral suspension be used in order to prevent clogging of the column.
- It is important to work quickly during this procedure.

### 1G. Cell Lysate Preparation from Viral Suspension

- a. Transfer up to 100  $\mu\text{L}$  of viral suspension to an RNase-free microcentrifuge tube (not provided).
- b. Add 350  $\mu\text{L}$  of **Lysis Solution**. Lyse viral cells by vortexing for 15 seconds. Ensure that mixture becomes transparent before proceeding to the next step.
- c. Add 200  $\mu\text{L}$  of 95 – 100% ethanol (provided by the user) to the lysate. Mix by vortexing for 10 seconds. **Proceed to Step 2 (page 9).**

## Section 2. Total RNA Purification from All Types of Lysate

**Note:** The remaining steps of the procedure for the purification of total RNA are the same from this point forward for all the different types of lysate.

### 2. Binding RNA to Column

- a. Assemble a column with one of the provided collection tubes
- b. Apply up to 600  $\mu\text{L}$  of the lysate with the ethanol (from **Step 1**) onto the column and centrifuge for 1 minute.

**Note:** Ensure the entire lysate volume has passed through into the collection tube by inspecting the column. If the entire lysate volume has not passed, spin for an additional minute.

- c. Discard the flowthrough. Reassemble the spin column with its collection tube.
- d. Depending on your lysate volume, repeat Step **2b** and **2c** as necessary.

### Optional Step:

Norgen's Total RNA Purification Kit isolates total RNA with minimal amounts of genomic DNA contamination. However, an optional **On-Column DNA Removal Protocol** is provided in Appendix A for maximum removal of residual DNA that may affect sensitive downstream applications. It is recommended that Norgen's RNase-Free DNase I Kit (Product # 25710) be used for this step. This step should be performed at this point in the protocol

### 3. Column Wash

- a. Apply 400  $\mu\text{L}$  of **Wash Solution** to the column and centrifuge for 1 minute.

**Note:** Ensure the entire wash solution has passed through into the collection tube by inspecting the column. If the entire wash volume has not passed, spin for an additional minute.

- b. Discard the flowthrough and reassemble the spin column with its collection tube.
- c. Repeat steps **3a** and **3b** to wash column a second time.
- d. Wash column a third time by adding another 400  $\mu\text{L}$  of **Wash Solution** and centrifuging for 1 minute.
- e. Discard the flowthrough and reassemble the spin column with its collection tube.
- f. Spin the column for 2 minutes in order to thoroughly dry the resin. Discard the collection tube.

### 4. RNA Elution

- a. Place the column into a fresh 1.7 mL Elution tube provided with the kit.
- b. Add 50  $\mu\text{L}$  of **Elution Solution** to the column.
- c. Centrifuge for 2 minutes at **200 x g (~2,000 RPM)**, followed by 1 minute at **14,000 x g (~14,000 RPM)** Note the volume eluted from the column. If the entire 50  $\mu\text{L}$  has not been eluted, spin the column at 14,000 x g (~14,000 RPM) for 1 additional minute.

**Note:** For maximum RNA recovery, it is recommended that a second elution be performed into a separate microcentrifuge tube (Repeat **Steps 4b** and **4c**).

### 5. Storage of RNA

The purified RNA sample may be stored at  $-20^{\circ}\text{C}$  for a few days. It is recommended that samples be placed at  $-70^{\circ}\text{C}$  for long term storage.

## Appendix A

### Protocol for Optional On-Column DNA Removal

Norgen's Total RNA Purification Kit isolates total RNA with minimal amounts of genomic DNA contamination. However, an optional protocol is provided below for maximum removal of residual DNA that may affect sensitive downstream applications. It is recommended that Norgen's RNase-Free DNase I Kit (Product # 25710) be used for this step.

1. For every on-column reaction to be performed, prepare a mix of 15  $\mu\text{L}$  of **DNase I** and 100  $\mu\text{L}$  of **Enzyme Incubation Buffer** using Norgen's RNase-Free DNase I Kit (Product # 25710). Mix gently by inverting the tube a few times. **DO NOT VORTEX.**

**Note:** If using an alternative DNase I, prepare a working stock of 0.25 Kunitz unit/ $\mu\text{L}$  RNase-free DNase I solution according to the manufacturer's instructions. A 100  $\mu\text{L}$  aliquot is required for each column to be treated.

2. Perform the appropriate Total RNA Isolation Procedure for your starting material up to and including **“Binding to Column”** (Steps 1 and 2 of all protocols)
3. Apply 400  $\mu\text{L}$  of **Wash Solution** to the column and centrifuge for 2 minute. Discard the flowthrough. Reassemble the spin column with its collection tube.
4. Apply 100  $\mu\text{L}$  of the RNase-free DNase I solution prepared in Step 1 to the column and centrifuge at 14, 000 x g (~14 000 RPM) for 1 minute.
 

**Note:** Ensure that the entire DNase I solution passes through the column. If needed, spin at 14, 000 x g (~14 000 RPM) for an additional minute.
5. After the centrifugation in Step 4, pipette the flowthrough that is present in the collection tube back onto the top of the column.
 

**Note:** Ensure Step 5 is performed in order to ensure maximum DNase activity and to obtain maximum yields of RNA, in particular for small RNA species.
6. Incubate the column assembly at 25 - 30°C for 15 minutes.
7. Without any further centrifugation, proceed directly to the second wash step in the **“Column Wash”** section (Step 3c).

## Appendix B

### Protocol for Total RNA Purification from Plasma or Serum

#### *Notes Prior to Use*

- Plasma or Serum of all human and animal subjects is considered potentially infectious. All necessary precautions recommended by the appropriate authorities in the country of use should be taken when working with plasma or serum.
- We recommend the use of this kit to isolate RNA from plasma or serum prepared by standard protocol from non-coagulating fresh blood using EDTA or sodium citrate as the anti-coagulant.
- Plasma prepared from fresh blood using heparin as an anti-coagulant is not suitable for use with this protocol. For heparin-prepared samples follow the protocol in section **1C, Lysate Preparation from Blood**.
- It is recommended that no more than 200  $\mu\text{L}$  of plasma or serum be used in order to prevent clogging of the column.
- Avoid multiple freeze-thaw cycle of the plasma or serum sample. Aliquot to the appropriate volume for usage prior to freezing.
- It is important to work quickly during this procedure.
- The yield of RNA from plasma and serum is highly variable. In general, the expected yield could vary from 1 to 100 ng per 100  $\mu\text{L}$  plasma or serum used. In addition, the expected A260:A280 ratio as well as the A260:A230 ratio will be lower (<1.80) than the normal acceptable range from other cells or tissues. Nonetheless, these isolated RNA could still be used effectively in different downstream applications such as RT-qPCR or microarrays.

#### 1. Cell Lysate Preparation from Plasma/Serum

- a. Transfer up to 200  $\mu\text{L}$  of plasma or serum to an RNase-free microcentrifuge tube (not provided).
- b. Add 300  $\mu\text{L}$  of Lysis Solution to every 100  $\mu\text{L}$  of plasma or serum. Mix by vortexing for 10 seconds.
- c. **Optional:** Add 0.7  $\mu\text{L}$  of 0.8  $\mu\text{g}/\mu\text{L}$  MS2 RNA per sample.

- Note:** The use of MS2 RNA could increase the consistency of downstream applications such as RT-PCR. However, the use of MS2 RNA is not recommended for applications involving global gene expression analysis such as microarrays or sequencing.
- d. Add 400  $\mu$ L of 95 – 100% ethanol (provided by the user) to every 400  $\mu$ L of the lysate (equivalent to every 100  $\mu$ L plasma or serum used). Mix by vortexing for 10 seconds. Proceed to Step 2 below.

## 2. Binding RNA to Column

- a. Assemble a column with one of the provided collection tubes
- b. Apply up to 600  $\mu$ L of the lysate with the ethanol (from Step 1) onto the column and centrifuge for 1 minute.

**Note:** Ensure the entire lysate volume has passed through into the collection tube by inspecting the column. If the entire lysate volume has not passed, spin for an additional minute.

- c. Discard the flowthrough. Reassemble the spin column with its collection tube.
- d. Depending on your lysate volume, repeat Step 2b and 2c as necessary.

### Optional Step:

Norgen's Total RNA Purification Kit isolates total RNA with minimal amounts of genomic DNA contamination. However, an optional On-Column DNA Removal Protocol is provided in Appendix A for maximum removal of residual DNA that may affect sensitive downstream applications. It is recommended that Norgen's RNase-Free DNase I Kit (Product # 25710) be used for this step. This step should be performed at this point in the protocol.

## 3. Column Wash

- a. Apply 400  $\mu$ L of 95 – 100% ethanol (provided by the user) to the column and centrifuge for 1 minute.

**Note:** Ensure the entire ethanol solution has passed through into the collection tube by inspecting the column. If the entire wash volume has not passed, spin for an additional minute.

- b. Discard the flowthrough and reassemble the spin column with its collection tube.
- c. Repeat steps 3a and 3b to wash column a second time.
- d. Wash column a third time by adding another 400  $\mu$ L of 95 – 100% ethanol (provided by the user) and centrifuging for 1 minute.
- e. Discard the flowthrough and reassemble the spin column with its collection tube.
- f. Spin the column for 2 minutes in order to thoroughly dry the resin. Discard the collection tube.

## 4. RNA Elution

- a. Place the column into a fresh 1.7 mL Elution tube provided with the kit.
- b. Add 50  $\mu$ L of Elution Solution to the column.
- c. Centrifuge for 2 minutes at 200 x g (~2,000 RPM), followed by 1 minute at 14,000 x g (~14,000 RPM) Note the volume eluted from the column. If the entire 50  $\mu$ L has not been eluted, spin the column at 14,000 x g (~14,000 RPM) for 1 additional minute.

**Note:** For maximum RNA recovery, it is recommended that a second elution be performed into a separate microcentrifuge tube (Repeat Steps 4b and 4c).

## 5. Storage of RNA

The purified RNA sample may be stored at  $-20^{\circ}\text{C}$  for a few days. It is recommended that samples be placed at  $-70^{\circ}\text{C}$  for long term storage.

## Troubleshooting Guide

Problem	Possible Cause	Solution and Explanation
Poor RNA Recovery	Incomplete lysis of cells or tissue	Ensure that the appropriate amount of Lysis Solution was used for the amount of cells or tissue.
	Column has become clogged	Do not exceed the recommended amounts of starting materials. The amount of starting material may need to be decreased if the column shows clogging below the recommended levels. See also "Clogged Column" below.
	An alternative elution solution was used	It is recommended that the Elution Solution supplied with this kit be used for maximum RNA recovery.
	Ethanol was not added to the lysate	Ensure that the appropriate amount of ethanol is added to the lysate before binding to the column.
	Ethanol was not added to the Wash Solution	Ensure that 50 mL of 96-100% ethanol is added to the supplied Wash Solution prior to use.
Poor RNA Recovery	Low RNA content in cells or tissues used	Different tissues and cells have different RNA contents, and thus the expected yield of RNA will vary greatly from these different sources. Please check literature to determine the expected RNA content of your starting material.
	Yeast: Lyticase was not added to the Resuspension Buffer	Ensure that the appropriate amount of lyticase is added when making the Resuspension Buffer.
	Bacteria and Yeast: All traces of media not removed	Ensure that all media is removed prior to the addition of the lysis solution through aspiration.
Clogged Column	Insufficient solubilization of cells or tissues	Ensure that the appropriate amount of lysis buffer was used for the amount of cells or tissue.
	Maximum number of cells or amount of tissue exceeds kit specifications	Refer to specifications to determine if amount of starting material falls within kit specifications
	High amounts of genomic DNA present in sample	The lysate may be passed through a 25 gauge needle attached to a syringe 5-10 times in order to shear the genomic DNA prior to loading onto the column.
	Centrifuge temperature too low	Ensure that the centrifuge remains at room temperature throughout the procedure. Temperatures below 15°C may cause precipitates to form that can cause the columns to clog.

Problem	Possible Cause	Solution and Explanation
RNA is Degraded	RNase contamination	RNases may be introduced during the use of the kit. Ensure proper procedures are followed when working with RNA. Please refer to "Working with RNA" at the beginning of this user guide.
	Procedure not performed quickly enough	In order to maintain the integrity of the RNA, it is important that the procedure be performed quickly. This is especially important for the Cell Lysate Preparation Step in the Tissue protocol, since the RNA in animal tissues is not protected after harvesting until it is disrupted and homogenized.
	Improper storage of the purified RNA	For short term storage RNA samples may be stored at $-20^{\circ}\text{C}$ for a few days. It is recommended that samples be stored at $-70^{\circ}\text{C}$ for longer term storage.
RNA is Degraded	Frozen tissues were allowed to thaw prior to RNA isolation	Do not allow frozen tissues to thaw prior to grinding with the mortar and pestle in order to ensure that the integrity of the RNA is not compromised.
	Starting material may have a high RNase content	For starting materials with high RNAase content, it is recommended that $\beta$ -mercaptoethanol be added to the Lysis Solution.
	Lysozyme or lyticase used may not be RNase-free	Ensure that the lysozyme and lyticase being used with this kit is RNase-free, in order to prevent possible problems with RNA degradation.
RNA does not perform well in downstream applications	RNA was not washed 3 times with the provided Wash Solution	Traces of salt from the binding step may remain in the sample if the column is not washed 3 times with Wash Solution. Salt may interfere with downstream applications, and thus must be washed from the column.
	Ethanol carryover	Ensure that the dry spin under the Column Wash procedure is performed, in order to remove traces of ethanol prior to elution. Ethanol is known to interfere with many downstream applications.
Genomic DNA contamination	Large amounts of starting material used	Perform RNase-free DNaseI digestion on the RNA sample after elution to remove genomic DNA contamination. It is recommended that Norgen's RNase-Free DNase I Kit (Product # 25710) be used for this step.

### Technical Support

Contact our Technical Support Team between the hours of 8:30 and 5:30 (Eastern Standard Time) at (905) 227-8848 or Toll Free at 1-866-667-4362.

Technical support can also be obtained from our website ([www.norgenbiotek.com](http://www.norgenbiotek.com)) or through email at [techsupport@norgenbiotek.com](mailto:techsupport@norgenbiotek.com).

### Product Use Restriction

Norgen's Total RNA Purification Kit Dx provides a rapid method for the isolation and purification of total RNA from tissue samples, blood, plasma, serum, bacteria, yeast, fungi and viruses for subsequent *in vitro* diagnostic use.

This kit is designed to be used with any downstream application employing enzymatic amplification or other enzymatic modifications of RNA followed by signal detection or amplification. Any diagnostic results generated using the RNA isolated with Norgen's Total RNA Purification Kit Dx in conjunction with an *in vitro* diagnostic assay should be interpreted with regard to other clinical or laboratory findings.

To minimize irregularities in diagnostic results, suitable controls for downstream applications should be used.

Norgen's Total RNA Purification Kit Dx is intended for use by professional users such as technicians, physicians and biologists experienced and trained in molecular biological techniques.


Norgen's Total RNA Purification Kit Dx does not provide a diagnostic result. It is the sole responsibility of the user to use and validate the kit in conjunction with a downstream *in vitro* diagnostic assay.

The respective user is liable for any and all damages resulting from application of Norgen's Total RNA Purification Kit Dx for use deviating from the intended use as specified in the user manual.

All products sold by Norgen Biotek are subjected to extensive quality control procedures and are warranted to perform as described when used correctly. Any problems should be reported immediately. The kit contents are for laboratory use only, and they must be stored in the laboratory and must not be used for purposes other than intended. The kit contents are unfit for consumption.

### Authorized Representative



 Norgen Biotek Corp.  
3430 Schmon Parkway, Thorold, ON Canada L2V 4Y6  
Phone: (905) 227-8848  
Fax: (905) 227-1061  
Toll Free in North America: 1-866-667-4362